

# Apollo

**TEGASUS**  
MEDICAL CONCEPTS INC.



## MEDICAL STORAGE CARTS

MEETS YOUR HYGIENE AND STORAGE DEMANDS

TAMBOUR DOOR

GLASS DOOR

SOLID DOOR

OPEN CART



# PEGASUS MEDICAL

**Pegasus Medical** is the partner of leading healthcare organizations and hospitals in 5 continents.

Our clients recommend our products and solutions since we greatly **simplify healthcare materials** management with flexible and scalable solutions.

We provide **significant cost savings**, and stand behind our methodology with **unmatched expertise and support**.

With over 20 years of experience with medication and hospital supply management, Pegasus Medical has developed our **Apollo carts** to provide an ideal solution for hospitals using ISO 60-40 trays and baskets.

# STORAGE CARTS

Pegasus Medical provides a unique, diversified storage cart product line. A combination of different enclosures, depths, heights and locking systems enables us to provide the optimum solution for any ward in your hospital or day care clinic.

Four solutions are available:

TAMBOUR DOOR

GLASS DOOR

SOLID DOOR

OPEN CART

Baskets are stored horizontally and can be tilted while pulled out





## OVERVIEW

Fully ISO (60X40cm) compatible

All surfaces are made of polymers and therefore resistant to moisture

ABS top cover, PP bottom floor, PVC exterior walls

The entire inside compartment of the cart is easy to clean

Self-drainage floor for quick and easy cleaning

5" (125mm) medical grade casters, non-marking rubber casters.

2x swiveling + 2x swiveling with brakes

Full 360° wrap-around bumper

Improved dense tray capacity due to the use of high capacity interior side tilted panels. Standard 183cm (72") cart accommodates up to 18 units of 5cm (2") deep trays / 11 pcs of 10cm (4") deep trays / 7 pcs of 20cm (8") deep trays

Extended caster pattern for increased stability

ABS side panels can be equipped with small wheels for smoother and easier movement of heavily loaded trays

Designated perforated rails accommodate the use of telescopic slides for support of heavy items. These can be used in combination with the side panels or separately

E Type



U Type



# TAMBOUR DOOR

TAMBOUR DOOR CARTS ARE EXTREMELY USEFUL  
IN SMALL ROOMS AND CORRIDORS WHERE  
THERE IS NOT ENOUGH SPACE FOR OPENING  
STANDARD DOORS



APOLLO E SINGLE CART  
APL-4.0E-TD



APOLLO UU DOUBLE CART  
APL-4.0UU-TD



## SMOOTH OPERATION


New counter-balanced system for smooth operation of roller shutter




## EASY CLEANING

Simple, tool-free removal of the roller shutter for cleaning purposes

## MULTIPLE DRAWER CONFIGURATIONS

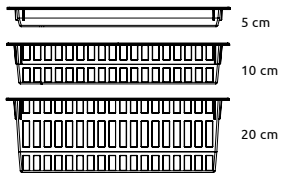

Modular drawer configurations, dividable baskets or solid trays

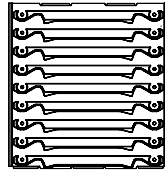


Tambour door Apollo carts are offered in different heights from 108cm (42.5") to 209cm (82.5")



### SECURITY & ERGONOMICS

A variety of locks e.g. code locks, key locks, network locks, offered on the top or middle part of the doors. Locks can be programmed remotely using the WIFI network enabling limited access to the staff

APOLLO UUU TRIPLE CART  
APL-4.0UUU-TD

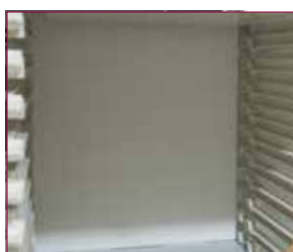

### PROTECTIVE WALL

Protective back wall prevents supplies from blocking the sliding tambour door



### ERGONOMICS

Mid-height handle eliminates the need to bend down to open the door

# GLASS DOOR

APOLLO CART LINE IS OFFERED WITH GLASS OR SOLID DOORS. GLASS DOORS ARE EASY TO OPERATE AND PROVIDE A CLEAR VIEW OF STORED PRODUCTS



APOLLO U SINGLE CART  
APL-4.0U-GD



APOLLO EE DOUBLE CART  
APL-2 .0EE-GD



## ERGONOMICS

Top to bottom handle for easier opening. Handle doesn't extend behind the footprint of the cart or protrude beyond




## DUST PROOF


Peripheral sealing to prevent dust accumulation, and reduce impact when the door is slammed



## MULTIPLE DRAWER CONFIGURATIONS

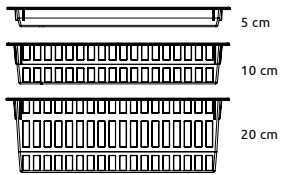

Modular drawer configurations, dividable baskets or solid trays

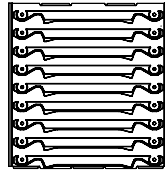

Glass door Apollo carts are offered in different heights from 108cm (42.5") to 209cm (82.5")



### Special ISO modules

A range of antimicrobial modules and high resistance modules for central sterilization departments is available upon request

ABS - Antimicrobial



PP - High resistance



ABS



PC - Transparent



APOLLO UUU TRIPLE CART  
APL-4.0-UUU-GD



### INFECTION CONTROL

All surfaces are designed for easy cleaning and draining. Carts can be offered in a machine washable version



### QUALITY GLASS

4mm thick tempered clear glass

# SOLID DOOR

SOLID DOORS ARE USED IN SPACES POPULATED WITH NON-MEDICAL STAFF, LIMIT PATIENTS' AND VISITORS' EXPOSURE TO THE MEDICAL STORED PRODUCTS



APOLLO E SINGLE CART  
APL-4.0E-SD



APOLLO UU DOUBLE CART  
APL-2.0UU-SD



## TOP WORKING SURFACE


Removable, washable ABS working surface with 3 sided gallery top and anti-spill lip on the front




## HPL SOLID DOORS

4mm HPL board  
High Pressure Laminate

## MULTIPLE DRAWER CONFIGURATIONS

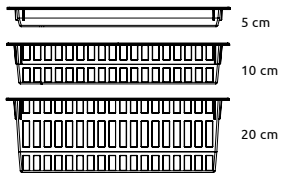

Modular drawer configurations, dividable baskets or solid trays

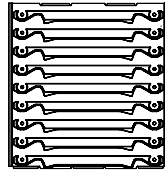

Solid door Apollo carts are offered in different heights from 108cm (42.5") to 209cm (82.5")



**APOLLO EE DOUBLE CART**  
APL-4.0EE-SD



### 270° OPENING

Unique hinge design enables 270° opening with a holding snap when attached to the cart side in an open position



### COLOR CODED LABELING HANDLE

Ergonomic tray handle improves visibility and increases nurses satisfaction

# OPEN CARTS

APOLLO CARTS WITHOUT DOORS ARE TYPICALLY USED IN ISOLATED ROOMS LIKE O.R., STERILE CORE WITH LOW DUST LEVELS OR LOCKED ROOMS WITH LIMITED ACCESS TO STAFF MEMBERS ONLY

The cassettes are according to ISO standard and fit Pegasus's different storage systems



APOLLO U SINGLE CART  
APL-4.0U-ND



APOLLO EE DOUBLE CART  
APL-4.0EE-ND




**QUALITY WHEELS**  
5" medical grade wheels produced by Tente




**BUMPERS**  
Both on the cart and on the room walls



## MULTIPLE DRAWER CONFIGURATIONS

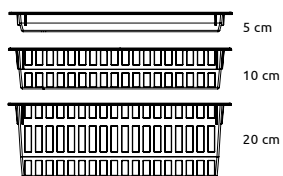

Modular drawer configurations, dividable baskets or solid trays

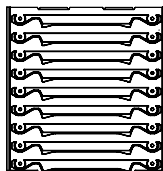

Open door Apollo carts are offered in different heights from 108cm (42.5") to 209cm (82.5")



ALL CARTS CAN BE FITTED IN ADVANCE WITH TELESCOPIC SLIDES OR ABS SIDE TILT PANELS OR A COMBINATION OF BOTH

\* For safety reason it is not allowed to load more than 10 kg in baskets placed on telescopic slides at upper levels (above 140cm)

APOLLO UUU TRIPLE CART  
APL-4.0UUU-ND



**DIN ACCESSORY RAIL**  
Allows easy installation of accessories




**HANDLES**  
Ergonomic lateral push handles

# CARTS SIZES


## E Type

E type carts are available as single, double or triple, with 4 height options. The depth is fixed at 59cm (23")

**Single cart E**  
75 cm (30")

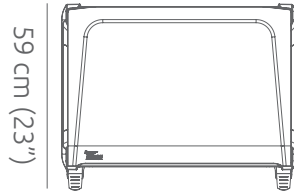


**Double cart EE**  
146 cm (58")

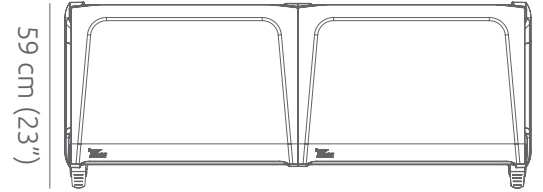


### 3 Panels cart

Height: 146 cm (57")  
APL-3.0E-TD/GD/SD/ND

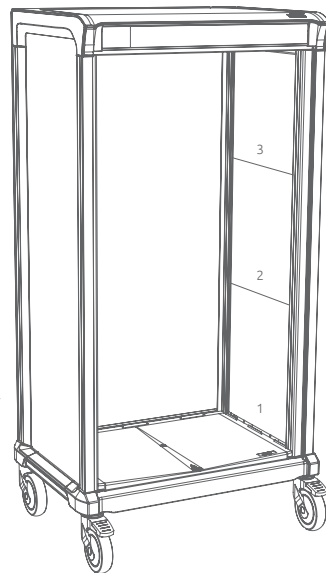


### 4 Panels cart

Height: 184 cm (72")  
APL-4.0E-TD/GD/SD/ND

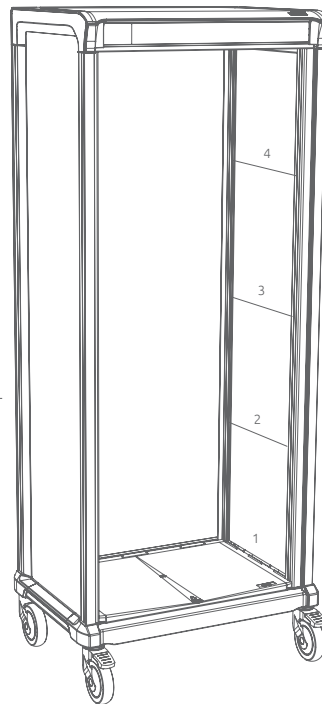


### 4.6 Panels cart

Height: 209 cm (82.5")  
APL-4.6E-TD/GD/SD/ND

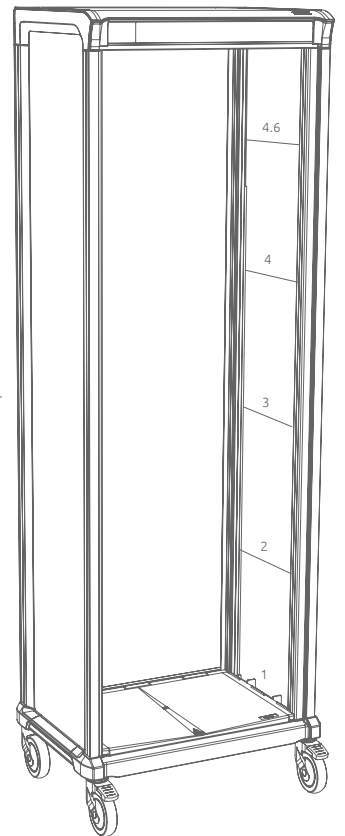

### 2 Panels cart

Height: 109 cm (43")  
APL-2.0E-TD/GD/SD/ND

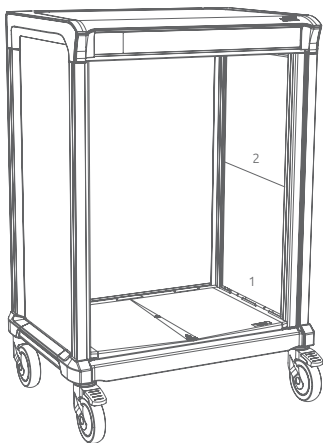

Trays or baskets maximum capacity for E and U type carts:


			2 PANELS	3 PANELS	4 PANELS	4.6 PANELS
20cm		TD	3	5	7	8
		GD/SD/ND	3	5	7	8
10cm		TD	5	8	11	13
		GD/SD/ND	6	9	12	14
5cm		TD	8	13	17	20
		GD/SD/ND	9	13	18	21

TD - Tambour Door    GD - Glass Door    SD - Solid Door    ND - No Door open cart


## U Type

U type carts are available as single, double or triple, with 4 height options. The depth is fixed at 79cm (31")


**Single cart U**  
55 cm (22")

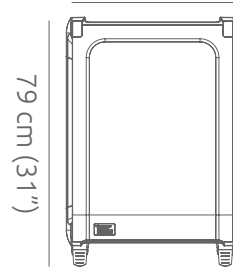


**Double cart UU**  
106 cm (42")

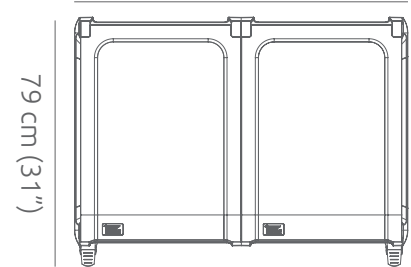

**4.6 Panels cart**  
Height: 209 cm (82.5")  
APL-4.6U-TD/GD/SD/ND

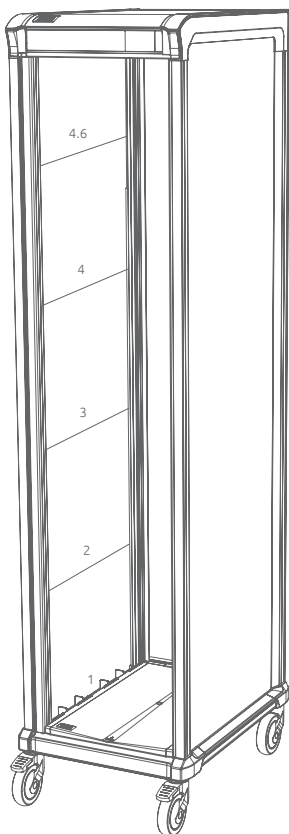


**4 Panels cart**  
Height: 184 cm (72")  
APL-4.0U-TD/GD/SD/ND

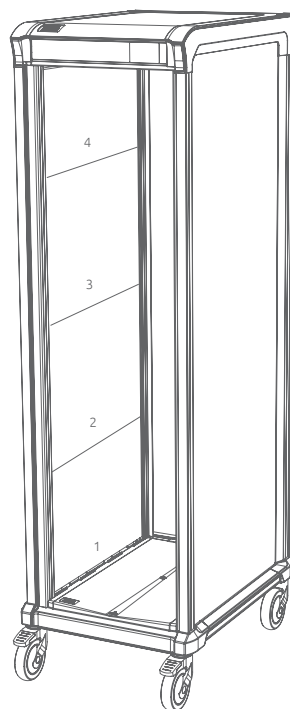

**3 Panels cart**  
Height: 146 cm (57")  
APL-3.0U-TD/GD/SD/ND



**2 Panels cart**  
Height: 109 cm (43")  
APL-2.0U-TD/GD/SD/ND

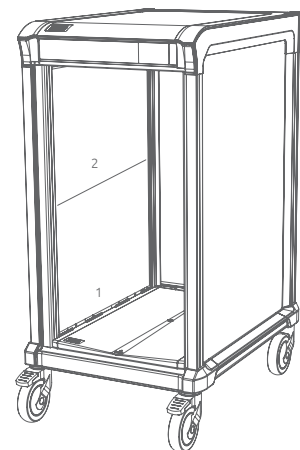



## CONTACT US FOR MORE DETAILS

+1 888-276-4750 | [info@pegasusmedical.net](mailto:info@pegasusmedical.net) | [www.pegasusmedical.net](http://www.pegasusmedical.net)

USA / CANADA  
905 E. Rose St Lakeland  
FL 33801 - USA

EUROPE / MIDDLE EAST / ASIA  
V. Croda Granda 13a 31027  
Spresiano (TV) - Italy